

-

Cardio Pulmonary Exercise Testing

Compact Metabolic cart for Accurate Pulmonary Gas Exchange & 12-lead Stress Testing ECG

Guark CPET

- Breath by Breath Gas Exchange data analysis (VO₂, VCO₂)
- Integrated 12-lead ECG for Stress Testing (option)
- Nutritional Assessment
- Fast response Paramagnetic O₂ Sensor
- Full Spirometry, Exercise SpO₂ monitor
- Mixing Chamber suitable for low and high ventilation ranges
- ► Windows Vista[™] Compatible



Applications

Designed to be the perfect tool for any kind of Cardio Pulmonary Exercise Testing (CPET), Quark CPET includes features for any discipline requiring metabolic assessment.

- Exercise Physiology, Education.
- Sports Science and Human Performance labs.
- Clinical Exercise Testing (Respiratory PFT Labs, Cardiology, Cardiac Rehab).
- Nutrition assessment.

Unsurpassed Accuracy

- Superior hardware specifications and quality design assure many years of accurate data.
- Quark CPET exceeds the accuracy standards of both the American Thoracic Society and the European Respiratory Society.
- Paramagnetic Oxygen and NDIR Carbon Dioxide analyzers.
- Bi-directional digital Turbine for Flow and Volume measurement.
- Accurate within a wide flow range (0-300 L/minute).
- Resistance to flow is less than 0.7 cmH,O/L/s@12 L/s.
- Small, Medium and Large sizes face masks for Adults.
- Comprehensive Calibration easy, quick and fully software assisted.
- Low Maintenance costs and Easy Servicing.

Cardio Pulmonary Exercise Testing


- ► True Windows[™] based program, offering, simplicity, versatility and familiar user interface.
- Display real-time data and graphs in either pre-defined or user formats.
- Real time O₂ and CO₂ waveforms.
- Exercise Flow-Volume loops.
- Ergometer control, via RS-232 interface, allows user easy protocol setup and dynamic changes.
- Advanced Data Elaboration.

- Automatic and manual detection of anaerobic threshold according to the Modified V-Slope method (Wasserman).
- Access data in a spreadsheet format for advanced data elaboration (filtering, smoothing, etc.).
- Standard and custom Exercise protocols design.
- O₂ Kinetics feature automatically provides O₂ debt, O₂ deficit and tau values during any constant stage.
- "Fitting" features on any plot (both linear and exponential).
- Indirect Cardiac Output by "Wassermann Algorithm".
- Export test in a single file (ASCII, Excel, file formats).
- Email test with MAPI compatible application (outlook, eudora etc.).
- Add user defined parameters and predicted equations with custom based formulas.
- Integrated 12-lead Stress Testing ECG.
- True diagnostic quality waveforms.
- Continuous 12 lead viewing.
- Zoom and freeze features.
- Current and reference ST analysis profiled for all 12 leads.
- ST depression and slope trends displayed during test.
- Averaged QRS complexes over reference ECG complex.
- Real time laser printout.



The CPET module add the possibility of performing pulmonary gas exchange analysis during exercise.



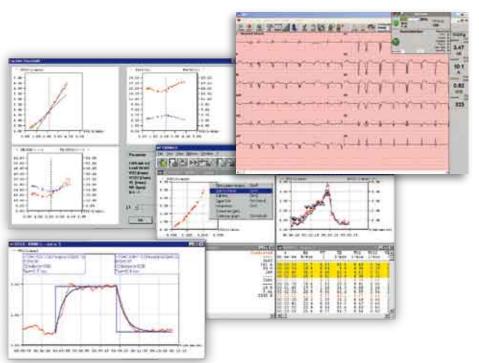

PC software Windows VISTA compatible offers familiar user interface and powerful data elaboration

Optional Hardware & Software

Mixing Chamber	Kit for low and high ventilation ranges 7 liters volume unit.	
Spirometry	Software and hardware to perform screening spi- rometry (FVC, MVV, bronchial challenge etc.).	
Pulse Oximetry	Monitor Oxygen saturation exercise using a pulse oximeter (SpO_2) provided with either finger, ear or reflectance probes.	
Canopy Kit	Perform REE measurements by means of canopy hood. Provided with specific flowmeter to in-crease accuracy at very low ventilations.	
High/Low FiO ₂ Kit	Kit of accessories for Gas Exchange measurements using hypoxic and hyperoxic gas mixtures.	
Ergometers	COSMED provides a wide selection of treadmills and bikes for both clinical and performance ap- plications. Refer to the Treadmill or Bike product brochures.	
Physiotrainer	Optional software to dynamically control the work load of any ergometer to target physiological "ef- fort-dependent" parameters such as VO ₂ , VO ₂ /Kg, HR and VE.	



12-lead stress test ECG available in both wireless and Patient-cable configurations .



Nutritional Assessment with canopy hood (option).

Nutritional Assessment


- REE, Substrate of Metabolism, NPRQ etc.
- User defined protocols consisting of two phases (1st phase discarded and 2nd phase data averaging).
- Post-analysis and custom selection of the "steady state" phase.
- Available w/ face mask and canopy hood (optional)

Easy Maintenance & Service

Quark CPET has been designed to minimize the need of ordinary and frequent service inspections. Built with the latest plug & play architecture service and support inquiries can now be handled not exclusively by either factory or specialized service sites.



Mlxing Chamber.

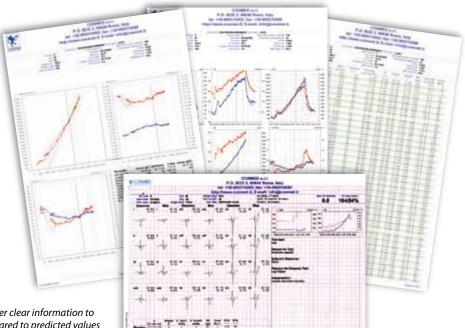

Powerful Reporting

- Prints any plot or data selection according to user defined criteria.
- Ready-to-print pre-defined reports.
- 9-plot Wasserman report can provide single page report with the 9 graphs and additional test results for an easy clinical interpretation.
- Summary Report, provides data for a simple and easy interpretation.
- ECG Report on standard or grid paper.

Explicative colour printout reports in different formats deliver clear information to user including: graphical test display, numerical data compared to predicted values and automatic interpretation of test results.

Technical Specification

Fest		Unit		
Pulmonary Gas Exchange (VO., VCO.)		•		
Anaerobic Threshold (AT)	2 2	•		
Indirect Cardiac Output (Was	sserman)	•		
Integrated Pulse Oximeter (SPO_)		0		
Integrated 12-lead ECG (GAS)	/EĈG)	0		
Physical Mixing Chamber		0		
Exercise Gas Exchange with E	levated FiO,	0		
HR Interface w/ external ECG (TTL)		•		
Resting Energy Expenditure (REE,RMR)		•		
Indirect Calorimetry (VO,, VCO., RQ)		•		
Indirect Calorimetry w/ Canopy Hood		0		
Indirect Calorimetry w/ Ventilator		0		
Forced/Slow Vital Capacity (SVC-FVC)		0		
Maximum Voluntary Ventilation (MVV)		0		
Broncho-Dilator test		0		
Broncho-Challenge Test		0		
nalyzers		Oxygen (O ₂)	Carbon dioxide (CO_2)	
Туре		Paramagnetic	NDIR	
Range		0-100%	0-10%	
Response time t90		<120 ms	<120 ms	
Accuracy		±0,03%	±0,03%	
Warm-up time		0 min	5 min	
Flowmeters		VO ₂ max	RMR (optional)	
Tuna		Bidirectional Digital Turbine	Bidirectional Digital Turbing	
Type		Ø 28 mm	Ø 18 mm	
Diameter		0.03-201/sec		
Flow Range		+ 2%	+ 2%	
Accuracy Resistance				
		<0.7 cm H2O s/L @ 12 L/s 5-300 I /m	< 0.7 cm H2O s/L @3 L/s	
Ventilation range		5-300 L/M	0-50 L/m	
nvironmental sensors	Temperature	Barometer	Humidity	
Range	0-50°C (32 - 122 F°)	400-800 mmHg	0-100%	
lardware				
Dimensions & weight		17 x 30 x 45 cm (6,2 x 12,9 x 16 in) / 6 Kg (13,2 lb)		
Dimensions & weight		17 X 30 X 13 CH1 (0)2 X 12)5 X 101	,,	
Dimensions & weight Interface ports		USB A-B, RS-232, HR-TTL, SPO,	.,,, =	

Gas Exchange Unit, HR belt, turbine flowmeter, PC software, 3 Adult face masks (S, M, L),5 mask disposable, 2 Adult Head caps, Calibration syringe (3 liters), Cables and probes, User manual

Available languages

Italian, English, German, Spanish, French, Portuguese.

Electrical requirements

Voltage	
Class	

100-240V ±10%; 50/60Hz 1 Type BF (EN 60601-1)

PC configuration required

Pentium or faster, Windows XP, VISTA 32, 128 Mb RAM or more, USB or RS 232, CD-Rom reader, 80 Mb on HD space available.

Safety & Quality Standards

Equipment complies with MDD (93/42 EEC) and FDA 510(k), EN 60601-1 (Safety) and EN 60601-1-2 (EMC)

COSMED and FITMATE names and logos are owned by COSMED and may be registered.





COSMED srl

Via dei Piani di Monte Savello 37 Pavona di Albano - Rome I - 00041 ITALY Phone +39 (06) 931-5492 Fax +39 (06) 931-4580 info@cosmed.com www.cosmed.com

COSMED USA Inc.

2211 N. Elston Avenue #305 Chicago, IL 60614 UNITED STATES Phone +1 (773) 645-8113 Fax +1 (773) 645-8116 info@cosmedusa.com www.cosmedusa.com

COSMED China Office

1st Floor, 215-1 QiYi Road Guangzhou 510030 P. R. of CHINA Phone +86 (20) 8332-4521 Fax +86 (20) 8332-0683 china@cosmed.it www.cosmed.com